


CAYA® Contoured Diaphragm

Intravaginal Barrier Device Product Description

Product Benefits:

- Easy insertion and removal
- Easy fit (one size fits most)
- Comfortable for both partners
- Eliminates latex-related odors

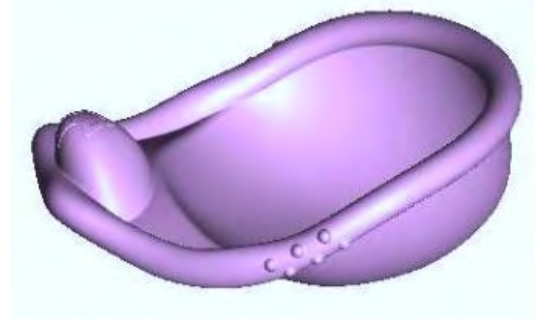

Product Specifications:

Elastic body:	50 shore-A durometer medical-grade silicone
Core:	Contoured, one-piece polymer spring
Overall length:	75 mm
Overall width:	67 mm
Membrane:	0.010" nominal thickness


Product Features:

Dynamic contours

- Folds in optimal insertion shape
- Embeds gently in supportive tissue

Gentle folding dynamics

- Soft, easy folding during insertion and removal
- Slides gently out of the vagina

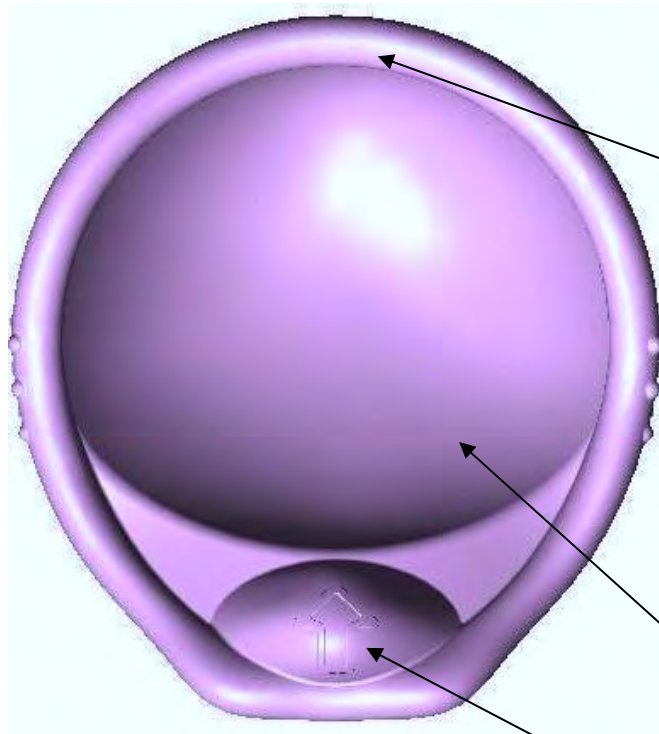

Firm insertion edge

- Folds compactly for easy insertion
- Stable shape when pushing into vagina
- Supports cervical cup in posterior fornix

Grip dimples (nubs)

- Enhanced grip when slippery
- Cue for bending locations

Cervical cup membrane

- Sized to fit a range of cervical positions
- Surrounds the cervix

Relief Arch

- Provides finger rest during insertion
- Extra clearance for partner

Finger-tip removal dome

- Enables easy hooking for removal (from top or bottom)
- Cue for insertion orientation